
THUASNE®

Lombatech

fr	Ceinture de soutien lombaire avec compression ajustable et ciblée.....	5
en	Lumbar support with adjustable and targeted compression	6
de	LWS-Stabilisierungssorthese mit anpassbarer und gezielter Kompression.....	7
nl	Lumbale rugbrace voor deonderrug met instelbare enverstevigde compressie	8
it	Cintura di sostegno lombare a compressione regolabile e mirata	9
es	Faja de sujeción lumbar con compresión ajustable y focalizada.....	10
pt	Cinta de suporte lombar com compressão ajustável e localizada	11
da	Lændestøttebælte med justerbar og målrettet kompression	12
fi	Mukautettava vahvistettu lannerangan tukivyö.....	13
sv	Ryggördel för ländryggen med justerbar och målinriktad kompression	14
el	Ζώνη οσφυϊκής υποστήριξης με ρυθμιζόμενη και στοχευμένη συμπίεση	15
cs	Bederní pás s cílenou a nastavitelnou kompresí	17
pl	Orteza lędźwiowa z regulowanym naciskiem ukierunkowanym	18
lv	Muguras jostas daļas atbalsta josta ar pielāgojamu un precīzu kompresiju	19
lt	Palaikomas reguliuojama tikslinių vietų suspaudimo juosmens diržas	20
et	Nimmepiirkonna reguleeritava ja kohandatava kompressiooniga tugikorsett.....	21
sl	Oporni križni pas z nastavljlivo in usmerjeno kompresijo.....	22
sk	Podporný driekový pás s nastaviteľnou a cielenou kompresiou	23
hu	Ágyéki fűző állítható és célzott kompresszióval	24
bg	Лумбален колан с регулируема и целенасочена компресия	25
ro	Centură de susținere lombară cu compresie reglabilă și direcționată	26
ru	Опорный поясничный корсет с регулируемой и направленной компрессией.....	27
hr	Lumbalni potporni pojas s podesivom i ciljanom kompresijom.....	29
zh	具有可调节向施压功能的腰部支撑带	30
ar	حزام لدعم أسفل الظهر بميزة الضغط القابل للتعديل والمستهدف	31

Instructions d'utilisation - Instructions for use - Gebrauchsanweisung - Gebruiksaanwijzingen - Istruzioni per l'uso - Instrucciones de uso - Instruções de utilização - Anvisninger vedrørende brug - Käyttöohjeet - Råd vid användning - Οδηγίες χρήσης - Pokyny k používání - Szczegóły użytkowania - Lietošanas instrukcija - Naudojimo instrukcijos - Kasutusjuhised - Navodila za uporabo - Pokyny na používanie - Használati utasítás - Инструкции за употреба - Instrucțiuni de utilizare - Инструкция по эксплуатации - Upute za uporabu - 使用指南 - إرشادات الاستعمال

	1	60 - 80 cm
	2	75 - 95 cm
	3	90 - 110 cm
	4	105 - 125 cm
	5	120 - 140 cm
	6	140 - 160 cm









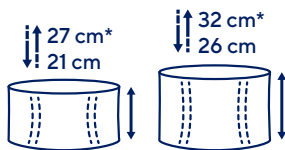
fr	Tour de taille	Hauteurs : 21 à 27* cm - 26 à 32* cm - *avec sangle additionnelle
en	Waist circumference	Heights: 21 to 27* cm - 26 to 32* cm - *with additional strap
de	Tailenumfang	Höhen: 21 bis 27* cm - 26 bis 32* cm - *mit Zusatzgurt
nl	Tailleomvang	Hoogtes: 21 tot 27* cm - 26 tot 32* cm - *met extra band
it	Circonferenza vita	Altezze: da 21 a 27* cm - da 26 a 32* cm - *con tirante aggiuntivo
es	Contorno de cintura	Alturas: 21 a 27* cm - 26 a 32* cm - *con correa adicional
pt	Perímetro de cintura	Alturas: 21 a 27* cm - 26 a 32* cm - *com fita adicional
da	Taljemål	Højde: 21 til 27* cm - 26 til 32* cm - *med ekstra strop
fi	Vyötäröympärys	Korkeudet: 21-27* cm - 26-32* cm - *lisähihnalla
sv	Midjemått	Höjd: 21 till 27* cm - 26 till 32* cm - *med extra rem
el	Περίμετρος μέσης	Ύψη: 21 έως 27* εκ. - 26 έως 32* εκ. - *με πρόσθετο ιμάντα
cs	Obvod pasu	Výšky: 21 až 27* cm - 26 až 32* cm - *spáslem dodatečného tahu
pl	Obwód w talii	Wysokości: 21 do 27* cm - 26 do 32* cm - *z dodatkowym paskiem
lv	Vidukļa apkārtmērs	Augstumi: 21 līdz 27* cm - 26 līdz 32* cm - *ar papildu siksnu
lt	Juostmens apimtis	Aukštis: 21-27* cm - 26-32* cm - *su papildomu diržu
et	Vööümbermõõt	Kõrgused: 21 cm kuni 27* cm - 26 cm kuni 32* cm - *lisarihmaga
sl	Obseg pasu	Višina: 21 do 27* cm - 26 do 32* cm - *dodatnim trakom
sk	Obvod pása	Výšky: 21 až 27* cm - 26 až 32* cm - *s prídavným popruhom
hu	Derékbőség	Magasság: 21 és 27* cm közötti - 26 és 32* cm közötti - *kiegészítő pánttal
bg	Обиколка на талията	Височини: от 21 до 27* cm - от 26 до 32* cm - *с допълнителна лента
ro	Circumferința taliei	Înălțimi: 21 până la 27* cm - 26 până la 32* cm - *cu fașă suplimentară
ru	Обхват талии	Высота: 21-27* см - 26-32* см - *с дополнительным ремнем
hr	Opseg struka	Visine: 21 dp 27* cm - 26 dp 32* cm - *s dodatnom postavljenom trakom
zh	腰围	高度：21至27*厘米 - 26至32*厘米 - *使用附加背带
ar	محيط الخصر	الارتفاعات: من 21 إلى 27* سم - من 26 إلى 32* سم - *بحزام إضافي

fr	Compression	Maintien	Éducation posturale
en	Compression	Support	Postural education
de	Kompression	Halt	Erlernen der Körperhaltung
nl	Compressie	Fixatie	Voorlichting over houding
it	Compressione	Sostegno	Educazione posturale
es	Compresión	Sujeción	Educación postural
pt	Compressão	Manutenção	Educação postural
da	Kompression	Støtte	Postural genoptræning
fi	Kompressio	Tuki	Asentokoulutus
sv	Kompression	Stöd	Postural träning
el	Συμπίεση	Συγκράτηση	Σωστή στάση σώματος
cs	Komprese	Stabilizace	Správné držení těla
pl	Kompresja	Podtrzymywanie	Edukacja dotycząca postawy
lv	Kompresija	Atbalsts	Stājas korekcija
lt	Kompresija	Palaikymas	Laikysenos taisymas
et	Kompressioon	Hooldus	Rühitreening
sl	Kompresija	Oprijem	Učenje telesne drže
sk	Kompresia	Podpora	Cvičenie posturálnych svalov
hu	Szorítás	Támasz	Tartásjavító oktatás
bg	Компресия	Поддръжка	Приучаване към правилна стойка
ro	Compresie	Susținere	Educație posturală
ru	Компрессия	Стабилизация	Исправление осанки
hr	Kompresija	Potpورا	Učenje o držanju
zh	压紧	稳定	姿势练习
ar	الضغط	الدعم	تدريب وضعي

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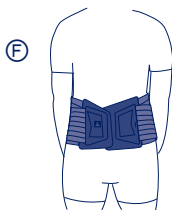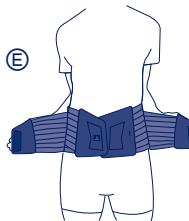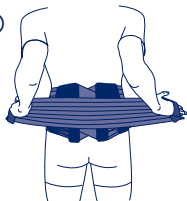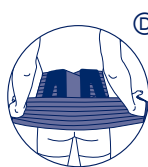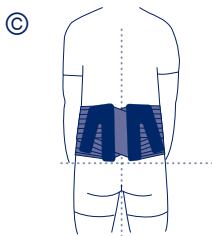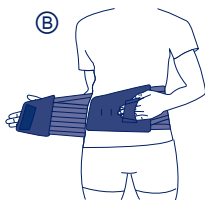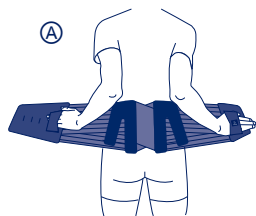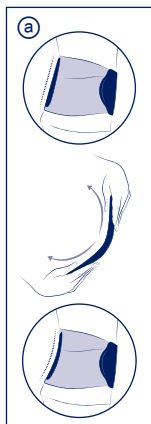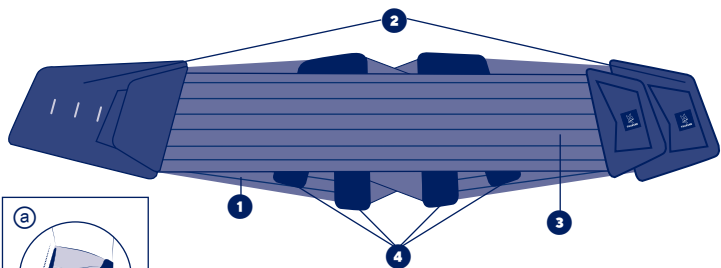Conserver cette notice.

en

LUMBAR SUPPORT WITH ADJUSTABLE AND TARGETED COMPRESSION

Description/Destination

The device is intended only for the treatment of the indications listed and for patients whose measurements correspond to the sizing table.

This device relieves lower back pain and stabilizes the spine, particularly when resuming activities.

Aerated and breathable Combitex fabric. ①

Abdominal panels with self-fastening surface and finger-loops. ②

Additional removable strap with finger-loops. ③

Comfortable dorsal stays. ④

Composition

Textile components: polyamide - polyester - elastane - polyurethane - acrylic resin.

Non-textile components: brushed tempered steel - polyoxymethylene - polypropylene - polyethylene.

Properties/Mode of action

The device relieves lower back pain by supporting the spine through an increase in intra-abdominal pressure. This reduces the load on the intervertebral discs.

Indications

Non-specific low back pain (acute, subacute and chronic phases).

Prevention of recurrent low back pain in combination with return to activity and/or therapeutic information. Spinal stenosis.

Contraindications

Do not use in the event of known allergy to any of the components.

Do not use on pregnant women.

Do not use in the event of a hiatus hernia.

Do not use in the event of bone cancer with spinal metastases.

Do not use in the event of circulatory, pulmonary or cardiovascular problems in patients for whom an increase in arterial blood pressure is not recommended.

Precautions

Verify the product's integrity before every use.

Do not use the device if it is damaged.

Wear the device above thin clothing.

Strictly comply with your healthcare professional's prescription and recommendations for use.

This product is intended for the treatment of a given condition. Its duration of use is limited to this treatment only.

Do not wear the device while sleeping.

In the event of discomfort, significant hindrance, pain, abnormal sensations, remove the device and consult a healthcare professional.

The use of a medical device by a child should be done under the supervision of an adult or a healthcare professional.

Measure the waist circumference of the patient at the mid-distance between the iliac crests and floating ribs.

Choose the appropriate size to fit the patient, referring to the size chart.

It is recommended that a healthcare professional supervises the first application.

For hygiene, security and performance reasons, do not re-use the product for another patient.

Do not wear the device in a medical imaging machine.

Undesirable side-effect

This device can cause skin reactions (redness, itching, irritations...).

Any serious incidents occurring related to the device should be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient is resident.

Instructions for use/Application

It is recommended that a healthcare professional adjust the stays to the shape of the patient's back before the belt is worn for the first time. ⑥

Position the belt around the waist, using the finger-loops located on the front panels. ⑦

The Thuasne® logo indicates which way up the belt should be put on.

Pull the finger-loops on the abdominal panels simultaneously and evenly.

Fasten the belt, right side over left side, adjusting to the desired compression. ⑧

The belt must be centred on the spine: the dorsal stays should be positioned on either side of the spine. ⑨

The bottom of the belt must stop at the top of the crease between the buttocks. ⑩

For a targeted side compression, use the additional strap. Position the additional strap on the painful area. ⑪

Pull the finger-loops on the additional strap simultaneously and evenly. ⑫

Attach the self-fastening parts of the additional strap onto the belt's abdominal panels or place them on top of each other. ⑬

The level of compression can be adjusted as required throughout the day.

Care

The product can be washed in accordance with the conditions shown on this instruction leaflet and on the label. Close the self-fastening tabs before washing. Machine washable at 30°C (delicate programme). It is not recommended to machine

wash the device more than 15 times (average: one wash per month). Do not dry clean. Do not use detergents, fabric softeners or aggressive products (products containing chlorine, etc.). Do not tumble-dry. Do not iron. Squeeze out any excess water. Dry away from any direct heat source (radiator, sun, etc.). Dry flat.

Storage

Store at room temperature, preferably in the original packaging.

Disposal

Dispose of in accordance with local regulations.

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Die Orthese muss zentral an der Wirbelsäule liegen: die Rücken-Stützstäbe müssen auf beiden Seiten der Wirbelsäule angebracht werden. ©
THUASNE UK Ltd
Thuasne SAS – SIREN/RCS Nanterre 542 091 186